# Angiotensin I (PRA)-RIA

*KIP5361*

# History

---

## Summary of change:

<b>Previous Version: 200116-1</b>	<b>Current Version: 200224-1</b>
No IVD logo	IVD logo added
	Addition of the following sentence at the end of the English IFU: "Other translations of this Instruction for Use can be downloaded from our website: <a href="https://www.diasource-diagnostics.com/">https://www.diasource-diagnostics.com/</a> "

Read entire protocol before use.

# Angiotensin I (PRA)-RIA

## I. INTENDED USE

Radioimmunoassay for the *in vitro* quantitative measurement of Angiotensin I in human EDTA plasma and the determination of Plasma Renin Activity (PRA) in human EDTA plasma.

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**For technical assistance or ordering information contact :**

**Tel : +32 (0) 10 84 99 00**

**Fax : +32 (0) 10 84 99 90**

## III. CLINICAL BACKGROUND

Renin is a proteolytic enzyme released by the juxtaglomerular cells of the kidney. The enzyme acts in general circulation to cleave its substrate, angiotensinogen, an alpha-2 globulin synthesized by the liver, in the decapeptide, angiotensin I.

Angiotensin I is a biologically inactive protein hormone and it is rapidly cleaved to a biologically active octapeptide angiotensin II, primarily in the lungs, by the activity of converting enzyme (ACE).

Angiotensin II has a very short half-life *in vivo* and is a very high potent vasopressor that plays an important role in hypertension and in blood pressure regulation. Angiotensin II plays also a role on aldosterone secretion by the adrenal gland.

The angiotensin II is rapidly degraded to inactive peptides fragments by enzymes present in plasma and tissues and for this reason it is difficult to measure the Angiotensin II levels in plasma.

As Angiotensin I represents directly the plasma renin activity, its measurement is used to determine the plasma renin activity (PRA) and to evaluate the renin-angiotensin system in disease states<sup>3,4</sup>

The PRA determination is useful for investigation of primary aldosteronism<sup>1,2</sup> (eg : adrenal carcinoma and adrenal cortical hyperplasia) and secondary aldosteronism<sup>5</sup> (renovascular disease, salt depletion, potassium loading, cardiac failure with ascites, pregnancy, Bartter syndrome)

The determination of PRA is also helpful to evaluate the prognosis and the appropriate therapy of patients with hypertension<sup>6</sup>.

Plasma renin activity (PRA) is quantified by measuring the rate of angiotensin I generation in plasma samples incubated at 37°C in presence of an enzyme inhibitor to avoid the conversion of angiotensin I to angiotensin II.

#### IV. PRINCIPLES OF THE METHOD

Plasma renin activity (PRA) is measured by the quantification of angiotensin I level in human plasma with a radioimmunoassay.

The method consists in 2 steps :

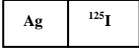
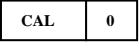
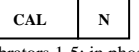
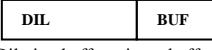
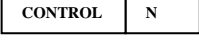
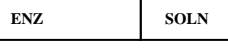
**Step 1:** Angiotensin I generation in plasma samples incubated at 37°C and at 4°C for the background value, in presence of an enzyme inhibitor to avoid angiotensin I enzymatic degradation.

**Step 2:** Radioimmunoassay :

A fixed amount of <sup>125</sup>I labelled angiotensin I competes with angiotensin I to be measured present in calibrators, controls and plasma samples for a fixed amount of antibody sites being immobilized to the wall of a polystyrene tube.

After 2 hours incubation at room temperature on a shaker an aspiration step terminates the competition reaction. The tubes are then washed with 2 ml of wash solution and aspirated again. A calibration curve is plotted and the angiotensin I concentrations of the samples are determined by dose interpolation from the calibration curve. Plasma renin activity of samples is then determined by subtracting the angiotensin I concentration found at 37°C by the angiotensin I concentration found at 4°C, considering the dilution factor used during enzymatic incubation and divided by the incubation time of enzymatic incubation .

#### V. REAGENTS PROVIDED

Reagents	96 Tests Kit	Colour Code	Reconstitution
 Tubes coated with anti angiotensin I	2 x 48	blue	Ready for use
 Ag <sup>125</sup> I <sup>125</sup> Iodine labelled angiotensin I HPLC grade	1 vial 150 kBq 10.5 ml	red	Ready for use
 CAL 0 Calibrator 0: in phosphate buffer with bovine serum albumin and Proclin	1 vial 2 ml	yellow	Ready for use
 CAL N Calibrators 1-5: in phosphate buffer with bovine serum albumin and Proclin (see exact values on vial labels)	5 vials 1 ml	yellow	Ready for use
 DIL BUF Dilution buffer: citrate buffer with Proclin and gentamycin	1 vial 12 ml	white	Ready for use
 WASH SOLN CONC Wash solution (TRIS-HCl)	1 vial 10 ml	brown	Dilute 70 x with distilled water (use a magnetic stirrer).
 CONTROL N Controls - N = 2 in phosphate buffer with bovine serum albumin and Proclin	2 vials 1 ml	silver	Ready for use
 ENZ SOLN Enzymatic solution: solution of enzymatic inhibitor in ethanol	1 vial 0.5 ml	white	Ready for use

**Note** :1.Use Calibrator 0 for dilution of samples with values above the highest calibrator before step 1.  
2.Calibrators were calibrated against NIBSC 86/536

#### VI. SUPPLIES NOT PROVIDED

The following material is required but not provided in the kit:

1. Distilled water
2. Pipettes for delivery of: 10 µl, 100 µl, 250 µl and 1 ml (the use of accurate pipettes with disposable plastic tips is recommended)

3. Glass tubes (12 x 75 mm) for enzymatic step (closed with a cap): 2 glass tubes needed for each plasma sample tested
4. Water bath at 37°C
5. Ice bath
6. Vortex mixer
7. Magnetic stirrer
8. Tube shaker (300 to 700 rpm)
9. 5 ml automatic syringe (Cornwall type) for washing
10. Aspiration system
11. Any gamma counter capable of measuring <sup>125</sup>I may be used (minimal yield 70%).

#### VII. REAGENT PREPARATION

**Working Wash solution:** Prepare an adequate volume of Working Wash solution by adding 69 volumes of distilled water to 1 volume of Wash solution (70x). Use a magnetic stirrer to homogenize. Discard unused Working Wash solution at the end of the day.

#### VIII. STORAGE AND EXPIRATION DATING OF REAGENTS

- Before opening, all kits components are stable until the expiry date, indicated on the label, if kept at 2 to 8°C.
- Freshly prepared working wash solution should be used on the same day.
- After its first use, tracer is stable, if kept in the original well-closed vial at 4°C until the tracer expiry date.
- After their first use Calibrators and controls are stable at 4°C until their expiry date .
- Alterations in physical appearance of kit reagents may indicate instability or deterioration.

#### IX. SPECIMEN COLLECTION AND PREPARATION

This kit is suitable only for EDTA plasma samples. Do not use heparin tubes. The following plasma collection procedure is recommended to prevent generation and degradation of angiotensin I.

Draw blood samples at room temperature. Do not chill EDTA collection tubes and don't store the collection tubes on ice.

Centrifuge the tube in a non-refrigerated centrifuge, and separate the plasma from the cells immediately after centrifugation, then aliquot and freeze immediately. Do not refrigerate. Cryoactivation, that is conversion of prorenin to renin, occurs when samples are chilled or exposed to temperatures of 6°C or lower for extended periods of time and when samples are chilled but remain in a liquid state. (i.e not frozen). Before analysis samples should be thawed rapidly to room temperature by standing them in lukewarm water.

Haemolyzed samples, icteric samples and lipemic samples must be discarded.

Avoid subsequent freeze-thaw cycles.

#### X. PROCEDURE

##### A. Handling notes

Do not use the kit or components beyond expiry date.  
Do not mix materials from different kit lots.

**Dilution buffer and enzymatic solution have to be stored at 2-8°C and pipetted directly after being taken out of the fridge (cf step B1 of the procedure)**

Bring all other reagents to room temperature prior to use.

Thoroughly mix all reagents and samples by gentle agitation or swirling.

Use a clean disposable pipette tip for addition of each different reagent and sample in order to avoid cross-contamination. High precision pipettes or automated pipetting equipment will improve the precision.

Prepare a calibration curve for each run, do not use data from previous runs.

##### B. Procedure

##### B1. Enzymatic Step (in glass tubes)

**Only for EDTA plasma samples, not for calibrators and controls**  
**This step needs a water bath at 37°C and an ice bath at 4°C**

1. For each EDTA plasma sample, label 2 glass tubes.
2. Place one tube in the ice bath.
3. Dispense 250µl of plasma sample, freshly taken from the fridge, in this tube.
4. Add 10 µl of enzymatic solution freshly taken from the fridge, into the tube.
5. Add 250 µl of dilution buffer freshly taken from the fridge, into the tube.
6. Vortex the tube.

- Transfer 250 µl of this mix in the second tube .
- Cover each tube with a cap.
- Place **immediately** one set of tubes into an ice bath at 4°C in a refrigerator, this is to determine the background value of angiotensin I.
- Place the second set of tubes into a water bath at 37°C, this is for determination of generated angiotensin I at 37°C.
- Incubate all tubes for 1 hour.
- After incubation, transfer immediately tubes incubated at 37°C into the ice bath to cool them and stop angiotensin I generation.

#### B2. Radioimmunoassay Step (in coated tubes)

- Label coated tubes in duplicate for each calibrator, control and treated plasma sample at 4°C and at 37°C. For the determination of total counts, label 2 normal tubes.
- Dispense 100 µl of calibrator or control or treated sample into the respective tubes
- Dispense 100 µl of <sup>125</sup>Iodine labelled angiotensin I into each tube, including the uncoated tubes for total counts .
- Incubate for 2 hours at room temperature (24 ± 4°C) on a tube shaker (300 to 700 rpm)
- Aspirate the content of each tube (except total counts). Be sure that the tip of the aspirator reaches the bottom of the coated tube in order to remove all the liquid.
- Wash tubes with 2 ml Working Wash solution (except total counts) and aspirate. Avoid foaming during the addition of the Working Wash solution.
- Let the tubes stand upright for two minutes and aspirate the remaining drop of liquid.
- Count tubes in a gamma counter for 60 seconds.

### XI. CALCULATION OF RESULTS

#### A. Calculation of angiotensin I concentrations

- Calculate the mean of duplicate determinations.
- Calculate the bound radioactivity as a percentage of the binding determined at the zero calibrator point (0) according to the following formula :

$$B/B_0 (\%) = \frac{\text{Counts (Calibrator or sample)}}{\text{Counts (Zero Calibrator)}} \times 100$$

- Plot the (B/B<sub>0</sub>(%)) values for each calibrator point as a function of angiotensin I concentration of each calibrator point. Reject obvious outliers.
- Computer assisted methods can be used to construct the calibration curve. If automatic result processing is used, a 4-parameter logistic function curve fitting is recommended.
- By interpolation of the sample (B/B<sub>0</sub> (%)) values, determine the angiotensin I concentrations of the samples from the calibration curve.
- For each assay, the percentage of total tracer bound in the absence of unlabelled angiotensin I (B<sub>0</sub>/T) must be checked.

#### B. Calculation of plasma renin activity (PRA)

Plasma renin activity (PRA) is expressed in ng/ml/hour of the in vitro generation of angiotensin I.

The background angiotensin I concentration measured on plasma samples at 4°C is subtracted from the angiotensin I concentration generated at 37°C, then multiplied by the initial sample dilution.

For the calculation of PRA by using the following equation :

$$PRA (\text{ng of angiotensin I /ml/ h}) = \frac{(\text{ng/ml at } 37^\circ\text{C} - \text{ng/ml at } 4^\circ\text{C}) \times 2.04}{\text{Enzymatic incubation time (h)}}$$

### XII. TYPICAL DATA

The following data are for illustration only and should never be used instead of the real time calibration curve.

Angiotensin I	cpm	B/Bo (%)
Total count	44052	
Calibrator		
0.0 ng/ml	24987	100.0
0.26 ng/ml	19438	77.8
0.75 ng/ml	11925	47.7
2.47 ng/ml	4420	17.7
9.11 ng/ml	1461	5.8
28.3 ng/ml	631	2.5

### XIII. PERFORMANCE AND LIMITATIONS

#### A. Detection Limit

The LOB (Limit of Blank) was calculated by measuring the blank 24 times and was calculated as the mean – 2 standard deviations of the distribution of these values. The LOB was calculated to be 0.06 ng/ml.

The LOD (Limit of Detection) was calculated as the LOB + 1.65 standard deviations of a low concentration sample. The LOD was calculated to be 0.14 ng/ml.

The LOQ (Limit of Quantification) was calculated by testing 5 samples of low values 10 times. The LOQ was calculated to be 0.28 ng/ml/h.

#### B. Specificity

The percentage of cross reaction was determined by testing plasma with spiked and unspiked crossreactants. The results are summarized in the following table:

Compound	Cross-Reactivity (%)
Angiotensin II	0.05
Angiotensin III	0.06
Angiotensinogen	ND
Tetradecapeptide	0.1

#### C. Precision

##### C1. Precision on Angiotensin I measurement without enzymatic step.

INTRA-ASSAY				INTER-ASSAY			
Sample	N	<X> ± SD (ng/ml)	C.V. (%)	Sample	N	<X> ± SD (ng/ml)	C.V. (%)
C1	20	1.46 ± 0.04	2.7	C1	14	1.38 ± 0.07	5.1
C2	20	4.37 ± 0.17	3.9	C2	14	4.61 ± 0.23	5.0

SD : Standard Deviation; CV: Coefficient of variation

##### C2. Precision on Plasma renin Activity calculated from tests on plasma measured in duplicates according the protocol of the kit.

INTRA-ASSAY				INTER-ASSAYS			
Sample	N	<X> ± SD (ng/ml/h)	C.V. (%)	Sample	N	<X> ± SD (ng/ml)	C.V. (%)
P1	10	1.99 ± 0.13	6.5	P1	6	0.34 ± 0.02	5.9
P2	10	10.35 ± 0.47	4.5	P14	6	0.42 ± 0.04	9.5

#### D. Accuracy

##### RECOVERY TEST

Plasma Sample	Added angiotensin I (ng/ml)	Measured angiotensin I (ng/ml)	Recovery (%)
1	8.28	9.09	109.7
	4.52	4.58	101.4
	2.75	2.35	85.5
2	8.28	7.94	95.9
	4.52	4.5	99.6
	2.75	2.55	92.8
3	7.74	6.17	79.7
	3.76	3.04	80.7
	1.64	1.59	96.8
4	7.74	6.75	87.2
	3.76	3.13	83.4
	1.64	1.62	99.3

## DILUTION TEST

Plasma samples were diluted with the zero calibrator, then tested following the procedure of the kit.

Plasma Sample	Sample dilution	Theoretical concent. (ng/ml)	Measured concent. (ng/ml)
1	1/1	15.36	15.36
	1/2	7.68	7.32
	1/4	3.84	3.32
	1/8	1.92	1.99
2	1/1	37.35	37.35
	1/2	18.67	18.13
	1/4	9.34	7.78
	1/8	4.67	4.14
	1/16	2.33	2.23
	1/32	1.17	1.35

### XIV. INTERNAL QUALITY CONTROL

If the results obtained for Control 1 and/or Control 2 are not within the range specified on the vial label, the results cannot be used unless a satisfactory explanation for the discrepancy has been given.

If desirable, each laboratory can make its own pools of control samples, which should be kept frozen in aliquots. Do not freeze-thaw more than twice.

Acceptance criteria for the difference between the duplo results of the samples should rely on Good Laboratory Practises.

### XV. EXPECTED VALUES

Each laboratory should establish its own range based on their local population. The PRA values here below are indicative only .

N	Patients	
19	Normal adult Supine	≤ 2.01 ng/ml/h
29	Normal adult Upright	≤ 4.95 ng/ml/h

### XVI. PRECAUTIONS AND WARNINGS

#### Safety

For *in vitro* diagnostic use only.

This kit contains <sup>125</sup>I (half-life: 60 days), emitting ionizing X (28 keV) and  $\gamma$  (35.5 keV) radiations.

This radioactive product can be transferred to and used only by authorized persons; purchase, storage, use and exchange of radioactive products are subject to the legislation of the end user's country. In no case the product must be administered to humans or animals.

All radioactive handling should be executed in a designated area, away from regular passage. A logbook for receipt and storage of radioactive materials must be kept in the lab. Laboratory equipment and glassware, which could be contaminated with radioactive substances, should be segregated to prevent cross contamination of different radioisotopes.

Any radioactive spills must be cleaned immediately in accordance with the radiation safety procedures. The radioactive waste must be disposed of following the local regulations and guidelines of the authorities holding jurisdiction over the laboratory. Adherence to the basic rules of radiation safety provides adequate protection.

The human blood components included in this kit have been tested by European approved and/or FDA approved methods and found negative for HbsAg, anti-HCV, anti-HIV-1 and 2. No known method can offer complete assurance that human blood derivatives will not transmit hepatitis, AIDS or other infections. Therefore, handling of reagents, serum or plasma specimens should be in accordance with local safety procedures.

All animal products and derivatives have been collected from healthy animals. Bovine components originate from countries where BSE has not been reported. Nevertheless, components containing animal substances should be treated as potentially infectious.

Avoid any skin contact with reagents (sodium azide as preservative). Azide in this kit may react with lead and copper in the plumbing and in this way form highly explosive metal azides. During the washing step, flush the drain with a large amount of water to prevent azide build-up.

Do not smoke, drink, eat or apply cosmetics in the working area. Do not pipette by mouth. Use protective clothing and disposable gloves.

### XVII. BIBLIOGRAPHY

- YOUNG WF Jr (1999)  
**Primary aldosteronism: A common and curable form of hypertension.**  
Cardiol Rev ; 7:207-214.
- YOUNG WF Jr (1997)  
**Pheochromocytoma and primary aldosteronism :diagnostic approaches.**  
Endocrinol Metab Clin North Am.,26:801-827.
- SEALEY J.E. (1975)  
**Radioimmunoassay of plasma renin activity.**  
Semin Nucl Med.,5(2):189-202.
- SPARANO F. (1978)  
**Behaviour of plasma renin activity and aldosterone during the first 72 h of life.**  
Clinical Endocrinology ,8, 207-211
- CUGINI P. (1992)  
**Secondary aldosteronism documented by plasma renin and aldosterone circadian rhythm in subjects with kidney or heart transplantation.**  
Ren Fail.,14(1):69-76
- W.GORDON WALKER (1976)  
**Relation between Plasma Renin Activity , Angiotensin, and Aldosterone and Blood Pressure in Mild Untreated Hypertension.**  
Circulation Research ,38(6),470-476

**XVIII. SUMMARY OF THE PROTOCOL**

	TOTAL COUNTS µl	CALIBRATORS CONTROLS µl	SAMPLE(S) µl
<b>ENZYMATIC STEP</b> (in one glass tube)			
Calibrators / controls	-	-	-
Samples	-	-	250
Enzymatic solution	-	-	10
Dilution buffer	-	-	250
Vortex tubes.  Transfer 250 µl of this mix in the second tube.  Cover each tube with a cap.  Place one tube into the ice bath in the refrigerator at 4°C and the other tube into the water bath at 37°C.  Incubate 1 hour.  Transfer tubes incubated at 37°C into the ice bath.			
<b>RADIOIMMUNOASSAY</b> (in coated tubes)			
Calibrators / controls	-	100	-
Samples from enzymatic step at 4°C and at 37°C	-	-	100
Tracer	100	100	100
Incubation	2 hours at 24±4°C on a shaker (300 to 700 rpm)		
Separation	-	Aspirate	
Working Wash solution	-	2.0 ml	
Separation	-	Aspirate	
Counting	Count tubes for 60 seconds		

Other translations of this Instruction for Use can be downloaded from our website: <https://www.diasource-diagnostics.com/>

DIAsource Catalogue Nr : KIP5361	Revision nr : 200224-1
-------------------------------------	---------------------------

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- YOUNG WF Jr (1999)  
**Primary aldosteronism: A common and curable form of hypertension.** *Cardiol Rev* ; 7:207-214.
- YOUNG WF Jr (1997)  
**Pheochromocytoma and primary aldosteronism :diagnostic approaches.** *Endocrinol Metab Clin North Am.*,26:801-827.
- SEALEY J.E. (1975)  
**Radioimmunoassay of plasma renin activity.** *Semin Nucl Med.*,5(2):189-202.
- SPARANO F. (1978)  
**Behaviour of plasma renin activity and aldosterone during the first 72 h of life.** *Clinical Endocrinology* ,8, 207-211
- CUGINI P. (1992)  
**Secondary aldosteronism documented by plasma renin and aldosterone circadian rhythm in subjects with kidney or heart transplantation.** *Ren Fail.*,14(1):69-76
- W.GORDON WALKER (1976)  
**Relation between Plasma Renin Activity , Angiotensin, and Aldosterone and Blood Pressure in Mild Untreated Hypertension.** *Circulation Research* ,38(6),470-476

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