

## Instruction for use

# EIA Mumps IgG

**REF** MuG096



Kit for professional use

**IVD** **CE**

 TestLine Clinical Diagnostics s.r.o.  
Křižíkova 68, 612 00 Brno, Czech Republic  
Tel.: +420 541 248 311  
FAX: +420 541 243 390  
E-mail: [info@testlinecd.com](mailto:info@testlinecd.com)  
[www.testlinecd.cz](http://www.testlinecd.cz)  
[www.testlinecd.com](http://www.testlinecd.com)



**CONTENT**

<b>1</b>	<b>Introduction .....</b>	<b>3</b>
<b>2</b>	<b>Test Principle.....</b>	<b>4</b>
<b>3</b>	<b>Materials Provided.....</b>	<b>5</b>
<b>4</b>	<b>Other Material Required for Manual Test Performance .....</b>	<b>6</b>
<b>5</b>	<b>Storage and Stability .....</b>	<b>6</b>
<b>6</b>	<b>Preparation of Reagents.....</b>	<b>6</b>
<b>7</b>	<b>Preparation of Samples .....</b>	<b>7</b>
<b>8</b>	<b>Assay Procedure.....</b>	<b>7</b>
<b>9</b>	<b>Working Schedule .....</b>	<b>9</b>
<b>10</b>	<b>Quality Control.....</b>	<b>10</b>
<b>11</b>	<b>Results Interpretation .....</b>	<b>11</b>
<b>12</b>	<b>Safety Precautions .....</b>	<b>12</b>
<b>13</b>	<b>Procedural Notes.....</b>	<b>12</b>
<b>14</b>	<b>IFU Symbols .....</b>	<b>14</b>

## Enzyme immunoassay for the detection of IgG antibodies to Mumps virus in human serum or plasma

---

### 1 Introduction

Mumps (parotitis epidemica) is an acute viral disease caused by an RNA virus from the *Paramyxoviridae* group.

Droplet-transmission is characteristic for the infections; saliva contaminated objects are less likely contributors to the spread of the disease. The disease occurs seasonally with the highest incidence in winter and spring, affecting mainly children aged 5-9 years and adolescents aged 15-19.

The incubation period of the disease is 2-3 weeks. The disease may occur asymptotically, with relatively mild non-specific symptoms (loss of appetite, elevated body temperature, headache), especially in children. A typical symptom of mumps is one-sided or bilateral painful swelling of parotid salivary glands, possibly also sublingual or submandibular salivary glands. Complications of the disease are more common in adults, the most severe being aseptic meningitis, orchitis, oophoritis and otitis.

The main means of mumps prevention is general vaccination of children using MMR vaccine, whereas the parotitic component of the vaccine induces the production of antibodies in about 90% of the vaccinated individuals.

The diagnosis of mumps is based on clinical picture and laboratory tests. Laboratory diagnostics is performed by direct detection of infectious agents or by the determination of specific antibodies using ELISA. Upon infection, specific antibodies are produced within 3-10 days after the onset of clinical signs. While IgM and IgA antibodies disappear after a few weeks or months, IgG antibodies persist for a long time, usually for life. Significant increases in IgG antibody levels occur after vaccination, although titres of these antibodies are generally lower than after natural infection and may not persist for life. The determination of IgG antibody levels also serves as a control of vaccination effect.

## 2 Test Principle

The kit is intended for detection of specific IgG antibodies in a sample by means of a sandwich type of the EIA method (i.e. a solid phase coated with specific antigen – antibody from the analysed sample – labelled antibody). The labelled antibody (conjugate) is an animal immunoglobulin fraction to human IgG conjugated with horseradish peroxidase. Peroxidase activity is determined in the test by a substrate containing TMB. Positivity is indicated when blue colour appears; after stopping solution has been added, blue changes to yellow. The yellow colour intensity is measured by a photometer at 450 nm, and it is proportional to the concentration of specific IgG antibodies in the sample.

### Antigen used

Purified and inactivated native antigen with high content of specific immunodominant epitopes.

### 3 Materials Provided

MICROPLATE	<b>Microtitre Plate</b>	1 pc
	coated with antigen, 12 x 8 wells in bag with desiccant	
CONTROL - CAL1	<b>Negative Control (Calibrator 1) 5 U/ml</b>	1 × 2 ml
	Solution containing no specific human antibodies, ready to use	
CUTOFF CAL2	<b>CUT-OFF (Calibrator 2) 20 U/ml</b>	1 × 3 ml
	Solution containing specific human antibodies in cut-off concentration, ready to use	
CONTROL + CAL3	<b>Positive Control (Calibrator 3) 80 U/ml</b>	1 × 2 ml
	Solution containing specific human antibodies, ready to use	
CAL4	<b>Calibrator 4 (320 U/ml)</b>	1 × 2 ml
	Solution containing specific human antibodies, ready to use	
CONJUGATE	<b>Conjugate</b>	1 × 15 ml
	Solution containing peroxidase labelled animal immunoglobulin to human IgG, ready to use	
DILUENT 2	<b>Sample Diluent 2</b>	1 × 105 ml
	Buffer with protein stabilisers, ready to use	
SUBSTRATE 2	<b>TMB-Complete 2</b>	1 × 15 ml
	Chromogenic substrate solution containing TMB/H <sub>2</sub> O <sub>2</sub> , ready to use	
WASH 20x	<b>Wash Solution</b>	1 × 75 ml
	20× concentrated buffer	
STOP	<b>Stop Solution</b>	1 × 15 ml
	Acid solution, ready to use	
	<b>Instructions for use</b>	1 pc

## 4 Other Material Required for Manual Test Performance

Single and multichannel pipettes

Disposable tips

Microplate washer

Timer

Incubator (37°C)

Microplate reader

## 5 Storage and Stability

Store the kit at +2°C to +8°C. Do not freeze. If the kit is stored as described, the labelled expiration date is valid. The expiration date is indicated on the package. The opened kit should be used within three months.

### Sample Preparation and Storage

The following human body liquids can be used for testing: serum and citrate plasma. Anticoagulants in the plasma (except for citrate) as well as bacterially contaminated, haemolytic or chylous samples can affect the test results.

Samples can be stored at +2°C to +8°C for one week. For a longer period, store samples at -20°C. Diluted samples should be used as soon as possible.

## 6 Preparation of Reagents

Dilute the Wash Solution 1:20 (1 part of solution and 19 parts of distilled water); e.g. 75 ml of the concentrated Wash Solution + 1425 ml of distilled water.

Salt crystals might develop in the bottle with the concentrated Wash Solution. Prior to use, it is necessary to dissolve the crystals by warming the bottle in a water bath. The diluted Wash Solution is stable at +2°C to +8°C for one week.

The Controls and the Calibrators are supplied ready to use, do not dilute further!

The Conjugate is supplied ready to use, do not dilute further!

TMB-Complete is a one-component chromogenic substrate solution ready to use, do not dilute further!

### Interchangeability of reagents

The Sample Diluent, TMB-Complete and the Avidity Solution are interchangeable in EIA kits of TestLine Clinical Diagnostics s.r.o., provided they have the identical numeric marking (e.g. Sample Diluent 2, Sample Diluent 3, etc.). The Stop Solution and the Wash Solution are universal in all kits.

## 7 Preparation of Samples

Mix gently the Sample Diluent prior to use.

### Dilution of sera and plasma samples

Dilute well mixed samples 1:101 with the Sample Diluent:

E.g.: 10 µl of sample + 1 ml of the Sample Diluent

Mix well.

## 8 Assay Procedure

Allow all reagents to come to room temperature and mix well. If you do not use a whole microplate, return unnecessary strips into the bag with desiccant. Seal the bag tightly and store at +2°C to +8°C. Keep dry!

1. Dispense the controls (calibrators) and the diluted samples according to the working schedule.

### Semiquantitative evaluation in Index of Positivity (IP)

- Leave A1 well empty (blank).
- Pipette 100 µl of the Negative Control (Calibrator 1) into 1 well.
- Pipette 100 µl of CUT-OFF (Calibrator 2) into 2 wells.
- Pipette 100 µl of the Positive Control (Calibrator 3) into 1 well.
- Pipette 100 µl of the diluted samples (see Chapter Preparation of Samples) into the other wells.

### Quantitative evaluation in Units U/ml

- Leave A1 well empty (blank).
  - Pipette 100 µl of the Negative Control (Calibrator 1) into 1 well.
  - Pipette 100 µl of CUT-OFF (Calibrator 2) into 2 wells.
  - Pipette 100 µl of the Positive Control (Calibrator 3) into 2 wells.
  - Pipette 100 µl of the Calibrator 4 into 2 wells.
  - Pipette 100 µl of the diluted samples (see Chapter Preparation of Samples) into the other wells.
2. Cover the microplate with the lid and incubate at 37°C for 30 minutes.
  3. Aspirate the content of the wells and wash 5× with the working strength Wash Solution. Fill the wells up to the edge. Finally, tap the inverted microplate thoroughly on an absorbent paper to remove solution remnants.
  4. Pipette 100 µl of the Conjugate into all wells except A1 well.
  5. Cover the microplate with the lid and incubate it at 37°C for 30 minutes.


6. Aspirate the content of the wells and wash 5× with the working strength Wash Solution. Fill the wells up to the edge. Finally, tap the inverted microplate thoroughly on an absorbent paper to remove solution remnants.
7. Pipette 100 µl of TMB-Complete into all wells. Avoid contamination – see Chapter Procedural Notes.
8. Cover the microplate with the lid and incubate at 37°C for 30 minutes. Keep out of light.
9. Stop the reaction by adding 100 µl of the Stop Solution in the same order and intervals as the substrate was added.
10. Read the colour intensity in wells against blank (A1 well) using photometer set to 450 nm. The absorbance should be read within 30 minutes after stopping the reaction.



### 9 Working Schedule

Semiquantitative evaluation

Index of Positivity (IP)

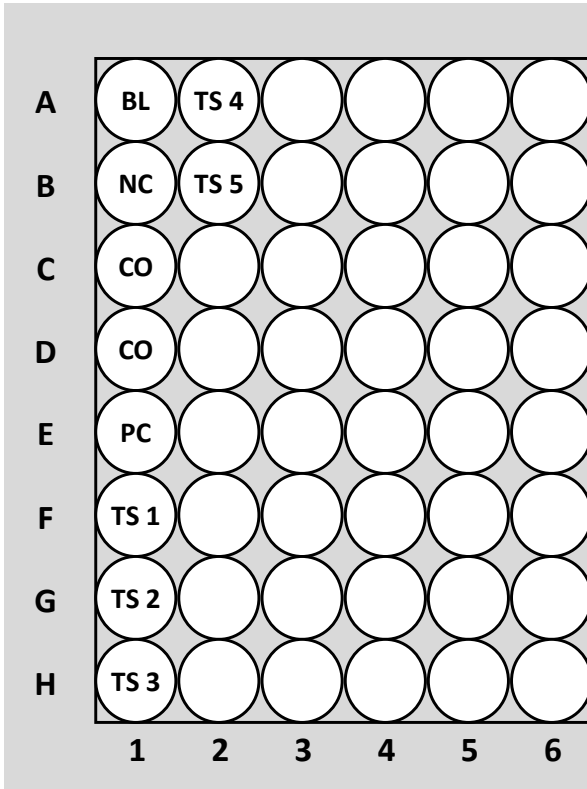

- BL      Blank (empty well)
- NC      100 µl    


CONTROL	-	CAL1
---------	---	------
- CO      100 µl    

CUTOFF	CAL2
--------	------
- PC      100 µl    

CONTROL	+	CAL3
---------	---	------
- TS 1-x   100 µl    diluted tested sample

Quantitative evaluation

Units U/ml

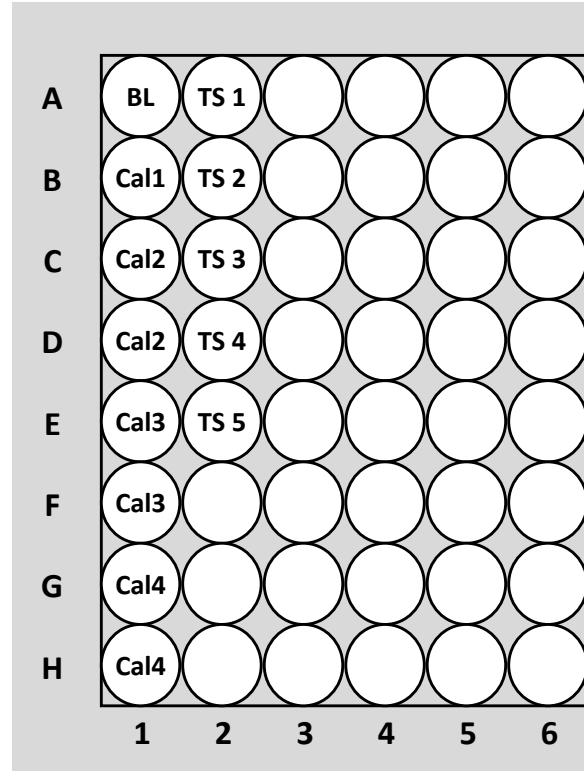

- BL      Blank (empty well)
- Cal1    100 µl    

CONTROL	-	CAL1
---------	---	------
- Cal2    100 µl    

CUTOFF	CAL2
--------	------
- Cal3    100 µl    

CONTROL	+	CAL3
---------	---	------
- Cal4    100 µl    

CAL4
------
- TS 1-x   100 µl    diluted tested sample

## 10 Quality Control

The test is valid if:

The absorbance of blank is lower than 0.150.

$$\text{BLANK} < 0.150$$

The absorbance of the Negative Control (Calibrator 1) is lower than half of the mean absorbance of CUT-OFF (Calibrator 2).

$$\boxed{\text{CONTROL}} - \boxed{\text{CAL1}} < 0.5 \times \boxed{\text{CUTOFF}} \boxed{\text{CAL2}}$$

The mean absorbance of CUT-OFF (Calibrator 2) is within a range of 0.150 – 0.900.

$$0.150 < \boxed{\text{CUTOFF}} \boxed{\text{CAL2}} < 0.900$$

The absorbance of the Positive Control (Calibrator 3) is 1.5-fold higher than the mean absorbance of CUT-OFF (Calibrator 2).

$$\boxed{\text{CONTROL}} + \boxed{\text{CAL3}} > 1.5 \times \boxed{\text{CUTOFF}} \boxed{\text{CAL2}}$$

The absorbance of the Calibrator 4 is higher than the absorbance of the Positive Control (Calibrator 3).

$$\boxed{\text{CAL4}} > \boxed{\text{CONTROL}} + \boxed{\text{CAL3}}$$

## 11 Results Interpretation

### Calculation of Index of Positivity (IP)

Divide the absorbance of a tested sample by the mean absorbance of CUT-OFF measured in the same test run:

$$IP = \frac{\text{Absorbance of sample}}{\text{Mean absorbance of CUT-OFF}}$$

Interpretation of the test results is described in Table 1.

**Table 1 Interpretation of test results**

Index of Positivity (IP)	Evaluation
lower than 0.9	negative
0.9 to 1.1	borderline
higher than 1.1	positive

Examination of borderline samples, i.e. samples with Index of Positivity from 0.9 to 1.1, should be repeated from a new sample collected after 2 to 6 weeks regarding to the disease specifics.

### Quantitative evaluation in Units (U/ml)

Construct a calibration curve by plotting the concentration (X) of the calibrators in U/ml against the corresponding absorbance (Y). Construct the calibration curve by single point cross connection. Read the values of antibody level (U/ml) in samples from the calibration curve. Interpretation of the quantitative test results is described in Table 2.

**Table 2 Quantitative interpretation in Units (U/ml)**

Antibody level (U/ml)	Evaluation
lower than 18	negative
18 to 22	borderline
higher than 22	positive

Examination of borderline samples should be repeated from a new sample collected after 2 to 6 weeks regarding to the disease specifics.

Serological finding can be interpreted only in the context of results of other laboratory tests and patient clinical picture.

## 12 Safety Precautions

The kit is intended for in vitro diagnostic use only.

The sera used for controls were tested and found to be negative for HIV 1 and HIV 2, HBsAg, HCV, TPHA. In spite of this fact, they still need to be handled as potentially infectious materials.

Some reagents contain sodium azide, which is a toxic compound. Avoid contact with skin.

The Stop Solution contains diluted acid solution. Avoid contact with eyes and skin.

It is necessary to observe the local safety rules and regulations.

### First aid

In case of contact with eyes, flush with copious amount of water and seek medical assistance. In case of contact with skin and clothing, remove all the contaminated clothes. Wash the skin with soap and plenty of running water. In case of contact with solutions containing plasma or clinical samples, disinfect the skin. In case of accidental ingestion, flush the mouth with drinking water and seek medical assistance.

### Remnants disposal

All the materials used for performing the test must be treated as potentially infectious due to the contact with biological materials. Therefore they need to be disposed together with biological waste.

### Expired kit disposal

Disassemble the kit and dispose the components as biological material. Discard the packaging material as required by local regulations.

## 13 Procedural Notes

In order to obtain reliable results, it is necessary to **strictly follow the Instructions for Use**. Always use clean preferably disposable tips and glassware.

**Microtitre Plate** – in order to prevent water condensation on the surface of the microplate, always allow the bag with the microplate to warm up to room temperature before opening.

**Wash Solution** – use high quality distilled water for preparing the working strength Wash Solution.

**Washing procedure** – keep to the prescribed number of wash cycles and fill the wells to the upper edge. The soak time (i.e. interval between two different wash cycles during which the wells stay filled up with the Wash Solution) should be approx. 30-60 seconds.

**TMB-Complete** – the vessel used for multichannel pipetting should not be used for other reagents. Do not return the surplus TMB-Complete from the pipetting vessel into the vial.

Non-reproducible results might be caused by improper methodology as following:

- insufficient mixing of reagents and samples before use
- improper replacement of vial caps
- using the same tip for pipetting different reagents
- reagent exposure to excessive temperature; bacterial or chemical contamination
- insufficient washing or filling of the wells (the wells should be filled to the upper edge), improper aspiration of Wash Solution remnants
- contamination of the well edges with Conjugate or samples
- using reagents from different kit lots
- contact of reagents with oxidants, heavy metals and their salts

The kit might be used for sequential examinations. When preparing working strength solutions, use only the amount of reagents needed for the analysis.


The kit might be used in all types of automatic EIA analysers.

If necessary, TestLine Clinical Diagnostics s.r.o. can offer a certified modification of the Instructions for Use for the specific type of analyser.

**The producer cannot guarantee that the kit will function properly if the assay procedure instructions are not strictly adhered to.**

## 14 IFU Symbols

---

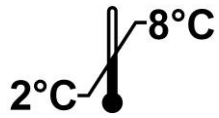

Temperature limitation

---



Keep dry

---



Expiry date

---



Lot number

---



Manufactured by

---



Consult instructions

---



Catalogue number

---



Number of tests

---














In vitro diagnostic medical device

---

---

Notes

**Summary of EIA Mumps IgG Protocol**

Step No.	Symbol	Test steps
1		Dilute samples serum/plasma 1:101 (10 µl + 1 ml)
2		Pipette Controls (Calibrators) and diluted samples – 100 µl Blank = empty well
3		Incubate at 37°C for 30 min
4		Aspirate and wash the wells 5×
5		Pipette Conjugate – 100 µl Blank = empty well
6		Incubate at 37°C for 30 min
7		Aspirate and wash the wells 5×
8		Pipette Substrate (TMB-Complete) – 100 µl Including blank
9		Incubate at 37°C for 30 min
10		Pipette Stop Solution – 100 µl Including blank
11		Read colour intensity at 450 nm